

Cardiac troponin T (cTnT)

The cardiac isoform of TnT is, similarly to cTnI, widely used as a marker of myocardial cell injury. cTnT has the same release kinetics into the bloodstream and the same sensitivity for minor myocardial injury as cTnI.

In human beings, cardiac troponin T is encoded by the TNNT2 gene. The major isoform found in normal adult human heart tissue (isoform 6 or TnT3) is 287 amino acids long with a calculated molecular weight of 34.6 kDa.

Reagents for immunoassay development

We provide MAbs that are suitable for the development of immunoassays for diagnostic purposes as well as several MAbs that are recommended for research use (see Figure 1). We also

provide polyclonal anti-cTnT antibodies as well as purified native and recombinant human cTnT and recombinant human slow and fast skeletal TnT proteins. The skeletal proteins are ideal for studying immunoassay cross-reactivity to these isoforms.

CLINICAL UTILITY

Early marker of acute myocardial infarction

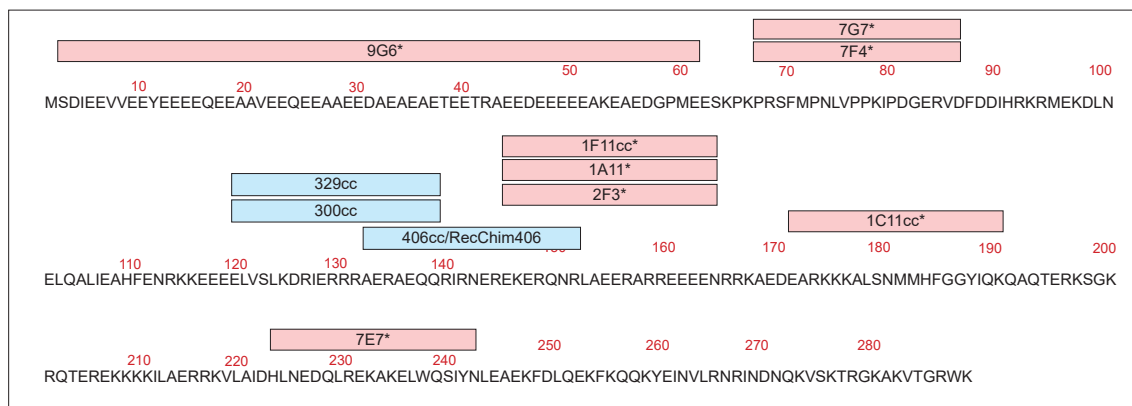


Figure 1. Epitope mapping of HyTest anti-cTnT monoclonal antibodies. We offer antibodies for the development of high-sensitivity cTnT assays (blue) as well as for research purposes (green, marked with *).

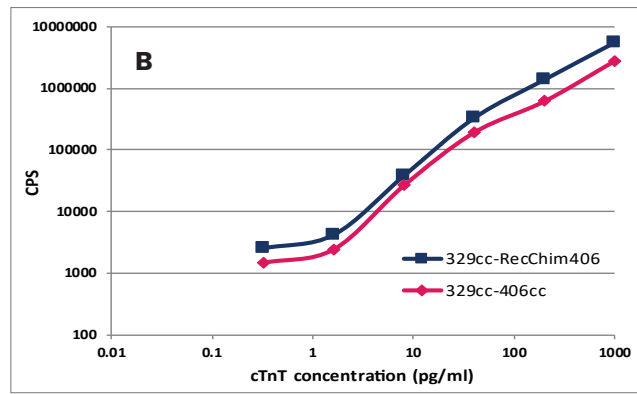
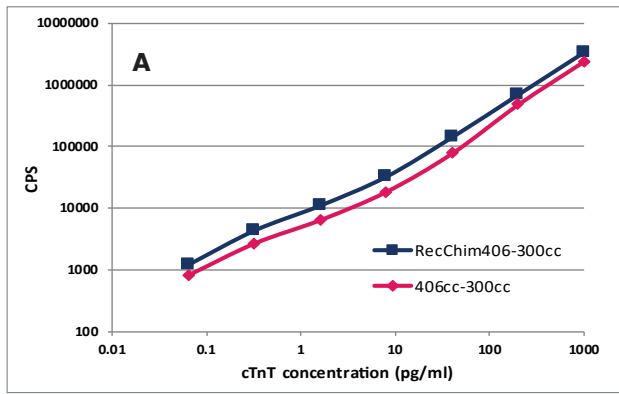


Figure 2. Comparing the performance of MAb 406cc and RecChim406 using a CLIA assay. In both assays 406cc or RecChim406 was used both as coating MABs (in pair with 300cc), and as biotin-labelled detection MABs (in pair with 329cc).

MONOCLONAL ANTIBODIES FOR HIGH-SENSITIVITY CTNT ASSAYS

We have developed three anti-cTnT MABs (300cc, 329cc and 406cc; Cat.# 4T19cc) that can be used for the development of an immunoassay with superior sensitivity (limit of detection better than 0.3 ng/l) and high specificity (no cross-reaction to cTnI or to skeletal isoforms of TnT up to 30 µg/l). MAB 406cc is also available as a recombinant chimeric construct in which the original wild type variable domains of the antibody and human IgG1 constant domains are combined (Cat.#RC4T19), MAB RecChim406). RecChim406 has the similar limit of detection

as the native antibody and like the latter could be used for the development of an immunoassay with superior sensitivity. In our preliminary tests, the RecChim406 is slightly more sensitive than the 406cc both as a capture and detection antibody (see Figure 2).

The ability of the antibody pairs 329cc-406cc and 406cc-300cc to recognize cTnT in the blood of AMI patients has been studied with over 80 serum and plasma samples. The antibody pairs demonstrate a good correlation with a commercially available hs-cTnT assay. Results of the analysis of 38 serum samples are provided in Figure 3.

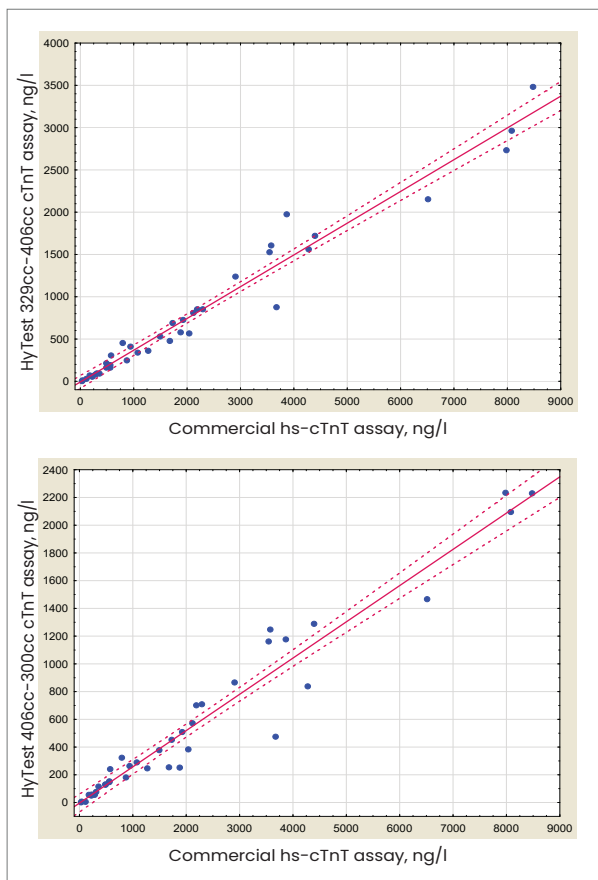


Figure 3. HyTest immunoassays show good correlation to a commercially available hs-cTnT assay. The concentration of cTnT in 38 serum samples obtained from AMI patients was determined by using two immunoassays that utilized HyTest antibodies (capture-detection pairs 329cc-406cc and 406cc-300cc) and a commercially available hs-cTnT assay.

Negligible cross-reactivity to skTnT

In high-sensitivity troponin assays, the specificity of the antibodies utilized is of utmost importance as even minor cross-reactivities could result in false positives.

We investigated the cross-reactivity of MABs 300cc, 329cc and 406cc to skeletal isoforms of troponin T. First, individual MABs were incubated with purified native cTnT and a mixture of recombinant slow and fast skTnT. All MABs recognized only cTnT (see Figure 4). Second, we tested the cross-reactivity of the two prototype assays to purified native skTnT, recombinant fast skTnT and recombinant slow skTnT. IGFBP-4 and MPO antigens were used as negative controls. Also in this case the cross-reactivity was well below 0.1% (see Figure 5).

Antibodies for research purposes

We offer several MABs that are recommended for research purposes. They also cross-react with cTnT proteins from different animal species (see Table 1).

Table 1. Cross-reactivity of anti-cTnT MABs with antigens from different animal species in Western blotting.

MAB	Human	Bovine	Porcine	Goat	Canine	Rabbit	Cat	Rat	Mouse	Fish
7F4	++	N/A	++	N/A	-	-	-	N/A	N/A	-
7G7	+	+	-	-	-	-	-	-	-	-
2F3	++	+	++	++	+	+	+	+	+	+
1A11	++	++	++	++	+	+	+	+	++	+
1F11	++	++	++	++	+	+	+	+	+	+

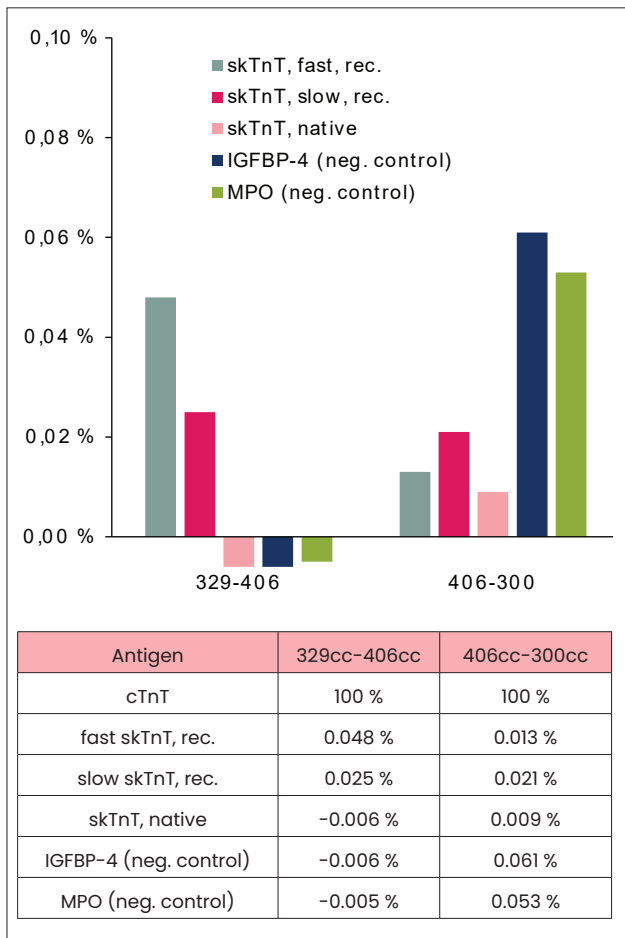
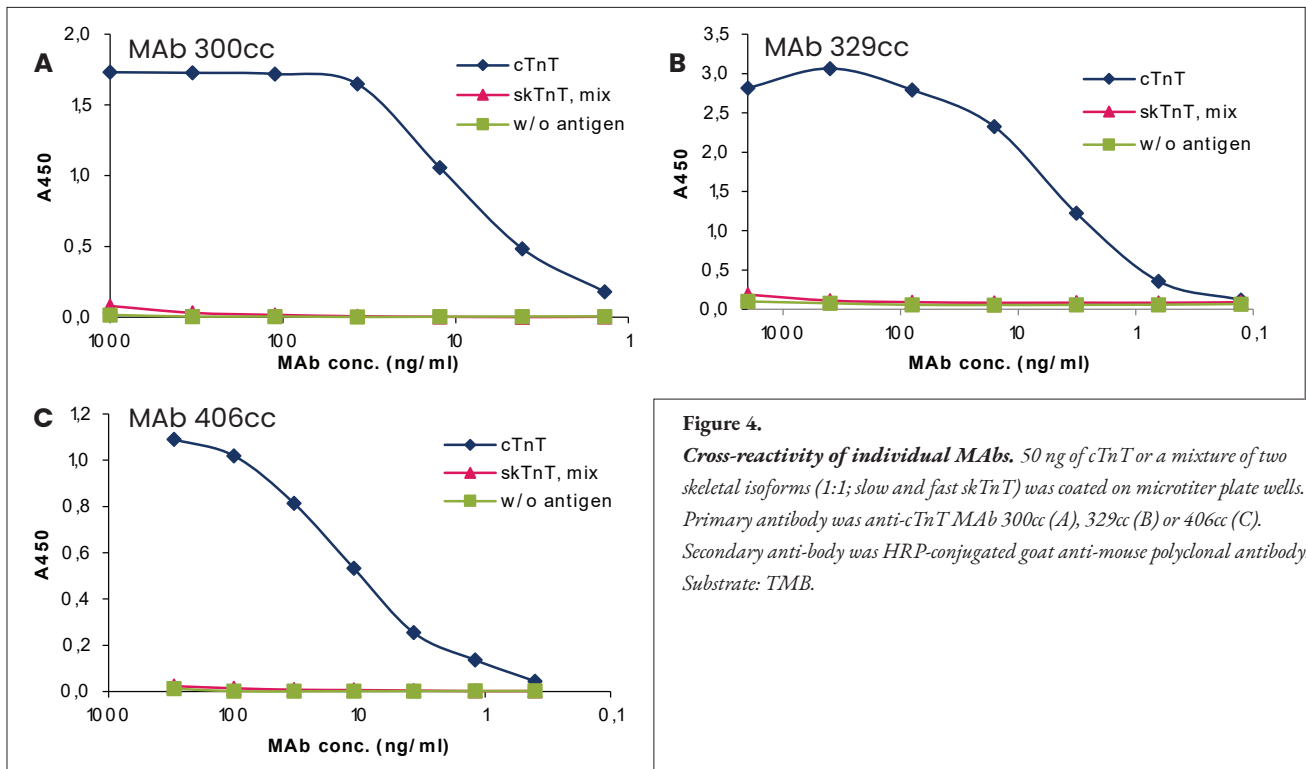


Figure 5.
Cross-reactivity of prototype assays 329cc-406cc and 406cc-300cc. Reactivity to various antigens (100 ng/ml) was investigated in sandwich immunoassays. Both assays demonstrated negligible cross-reactivity to all markers tested. In bs-cTnT immunoassays the cross-reactivity should be < 0.10 % and preferably < 0.05-0.03 %.

Chimeric antibodies prevent the HAMA effect

The performance of different combination of chimeric and original mouse antibodies have been tested using human anti-mouse antibody (HAMA) containing serum samples obtained from different individual HAMA-positive donors in order to verify the sensitivity of original mouse and chimeric antibodies to the HAMA effect. All combined mouse-and-chimeric pairs showed either no or just a negligible background signal with all tested HAMA samples (see Figure 6).

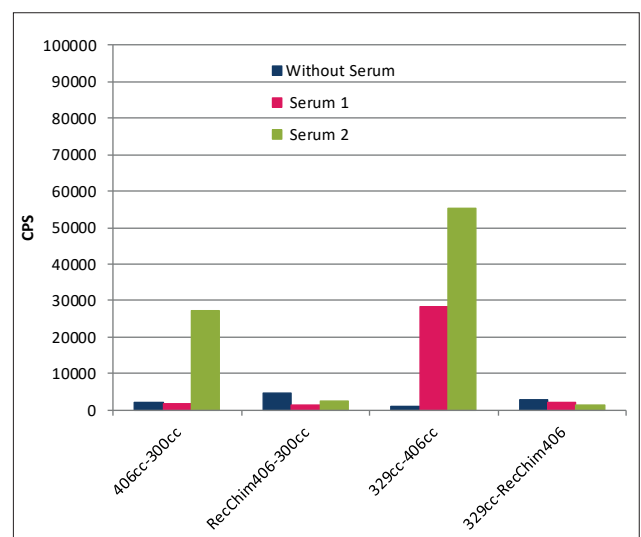


Figure 6.
Chimeric antibodies eliminate the HAMA effect. The performance of chimeric and native 406 in the presence of HAMA was tested in CLIA assay with two serum samples with the following HAMA concentrations: 807 ng/ml in sample 1; and 1,388 ng/ml in sample 2. As a control, buffer without serum was used.

PURIFIED ANTIGENS

Recombinant human slow and fast skTnT

The recombinant slow skeletal TnT (Cat.# 8RST2) and fast skeletal TnT (Cat.# 8RFT4) are ideal for studying immunoassay cross-reactivity to these isoforms.

Recombinant human cTnT

Isoform 6 (which is also known in the literature as TnT3) is the major isoform of troponin T that is presented in normal adult human heart tissue.

Our recombinant human cTnT (Cat.# 8RTT5) is produced in *E. coli* by expressing a gene encoding for the 288 amino acid long isoform 6 (TnT3) of cTnT. This isoform is the main isoform of cTnT in normal adult human heart tissue. The protein has an additional Met residue at its N-terminus.

ORDERING INFORMATION

MONOCLONAL ANTIBODIES

Product name	Cat. #	MAB	Subclass	Remarks
Troponin T cardiac	4T19	9G6	IgG1	EIA, WB, a.a.r. 2-61
		7F4	IgG2b	EIA, WB, a.a.r. 67-86
		7G7	IgG1	EIA, WB, a.a.r. 67-86
		2F3	IgG2b	EIA, WB, a.a.r. 145-164
		1A11	IgG2b	EIA, WB, a.a.r. 145-164
		7E7	IgG1	EIA, WB, a.a.r. 223-242
	4T19cc	300cc	IgG1	<i>In vitro</i> , EIA, a.a.r. 119-138
		329cc	IgG1	<i>In vitro</i> , EIA, a.a.r. 119-138
		406cc	IgG2a	<i>In vitro</i> , EIA, a.a.r. 132-151
		1F11cc	IgG2b	<i>In vitro</i> , EIA, WB, a.a.r. 145-164
		1C11cc	IgG1	<i>In vitro</i> , EIA, WB, a.a.r. 171-190
	RC4T19	RecChim406	IgG1	EIA, recombinant chimeric antibody

HUMAN ANTIGENS

Product name	Cat. #	Purity	Source
Troponin T cardiac, human, recombinant	8RTT5	>95%	Recombinant
Troponin T skeletal muscle, human	8T24	>95%	Human skeletal muscle
Troponin T fast skeletal, human, recombinant	8RFT4	>95%	Recombinant
Troponin T slow skeletal, human, recombinant	8RST2	>95%	Recombinant
Troponin complex (I-T-C), human	8T62	N/A	Human cardiac muscle
Troponin complex (I-T-C), artificial	8T62a	N/A	Human cardiac muscle